

Arterial Blood Samplers

XQZ-01 Preset Type



④ Rubber Block:

To seal the needle, preventing needle stick injury.



④ Luer Cap:

To seal the sample, ensuring sample quality.

④ Luer Lock Connection:

To prevent needle from falling off.

④ Vented Plug:

Double-layer porous structure enables rapid air venting while preventing blood leakage.


④ Preset System:

Automatic sampling of blood, simple and efficient.



④ Anticoagulant:

Ultra-precision sprayed solid heparin lithium, does not dilute blood samples. Balanced by calcium ions, avoids the combination of calcium ions and heparin in the blood, ensure the accuracy of the test results.

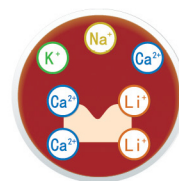


MOD NO.	Syringe Size	Needle Size (Gauge × Length, inch)	Blood Preset Volume	Heparin Lithium IU
XQZ-01 1mL 0.5mm×16mm	1 mL	25G × 5/8"	0.6mL	30IU
XQZ-01 1mL 0.5mm×25mm	1 mL	25G × 1"	0.6mL	30IU
XQZ-01 1mL 0.6mm×25mm	1 mL	23G × 1"	0.6mL	30IU
XQZ-01 1mL 0.7mm×25mm	1 mL	22G × 1"	0.6mL	30IU
XQZ-01 1mL 0.7mm×32mm	1 mL	22G × 1 1/4"	0.6mL	30IU
XQZ-01 3mL 0.5mm×16mm	3 mL	25G × 5/8"	1.6mL	80IU
XQZ-01 3mL 0.6mm×25mm	3 mL	23G × 1"	1.6mL	80IU
XQZ-01 3mL 0.7mm×25mm	3 mL	22G × 1"	1.6mL	80IU
XQZ-01 3mL 0.7mm×32mm	3 mL	22G × 1 1/4"	1.6mL	80IU

Shandong Ande Healthcare Apparatus Co., Ltd.

Arterial Blood Samplers

XQZ-02 Preset & Safety Type

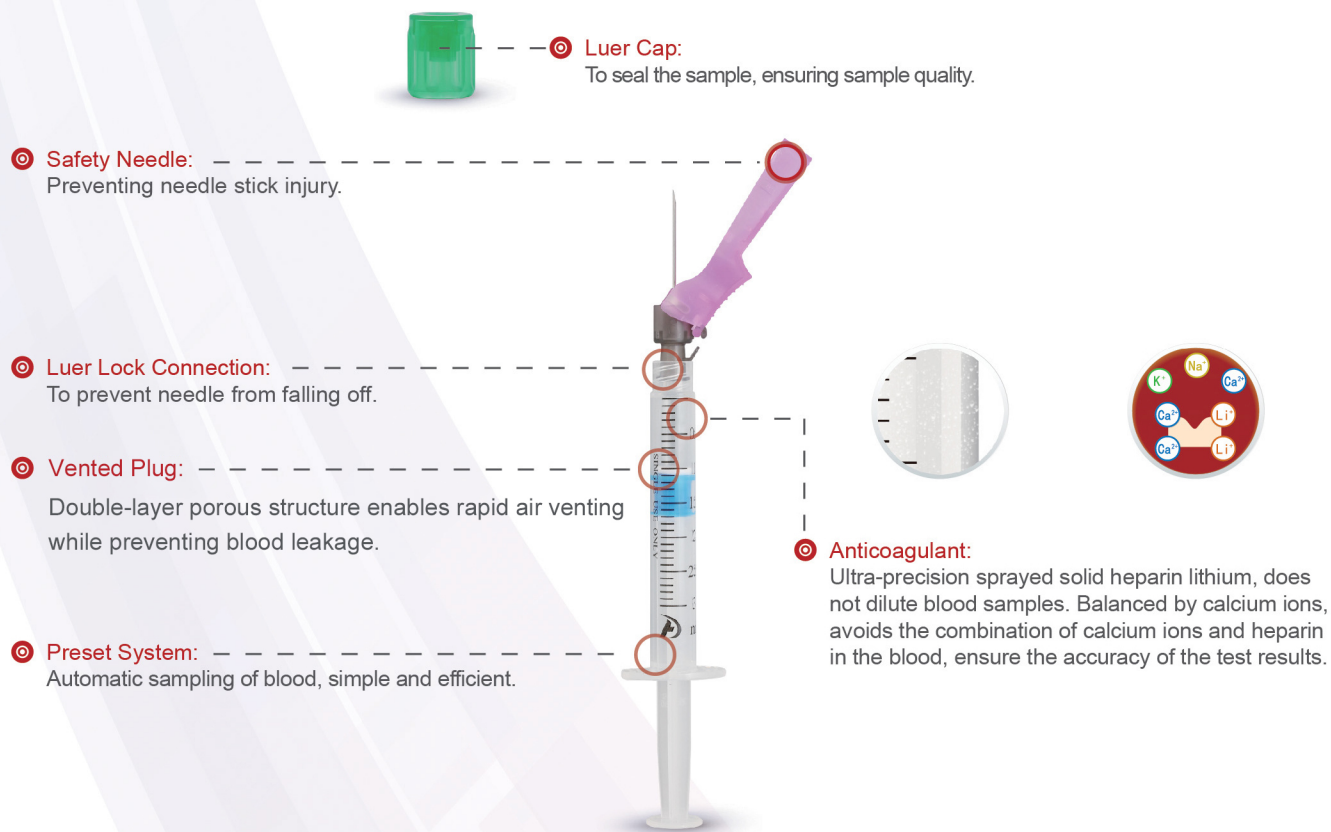

MOD NO.	Syringe Size	Needle Size (Gauge × Length, inch)	Blood Preset Volume	Heparin Lithium IU
XQZ-02 1mL 0.5mm×16mm	1 mL	25G × 5/8"	0.6mL	30IU
XQZ-02 1mL 0.5mm×25mm	1 mL	25G × 1"	0.6mL	30IU
XQZ-02 1mL 0.6mm×25mm	1 mL	23G × 1"	0.6mL	30IU
XQZ-02 1mL 0.6mm×32mm	1 mL	23G × 1 1/4"	0.6mL	30IU
XQZ-02 1mL 0.7mm×25mm	1 mL	22G × 1"	0.6mL	30IU
XQZ-02 1mL 0.7mm×32mm	1 mL	22G × 1 1/4"	0.6mL	30IU
XQZ-02 3mL 0.5mm×16mm	3 mL	25G × 5/8"	1.6mL	80IU
XQZ-02 3mL 0.6mm×25mm	3 mL	23G × 1"	1.6mL	80IU
XQZ-02 3mL 0.6mm×32mm	3 mL	23G × 1 1/4"	1.6mL	80IU
XQZ-02 3mL 0.7mm×25mm	3 mL	22G × 1"	1.6mL	80IU
XQZ-02 3mL 0.7mm×32mm	3 mL	22G × 1 1/4"	1.6mL	80IU

Shandong Ande Healthcare Apparatus Co., Ltd.

Arterial Blood Gas Syringes

XQZ-03 Without Needle



MOD NO.	Syringe Size	Fill Volume	Heparin Lithium IU
XQZ-03 1mL	1 mL	0.6mL	30IU
XQZ-03 3mL	3 mL	1.6mL	80IU

Shandong Ande Healthcare Apparatus Co., Ltd.