

PERIPHERALLY INSERTED CENTRAL VENOUS CATHETERS

PICC Contents

- Distal Valve PU Material
- Midline Catheter
- CT USE
- Neonate 1.9F



SPECIFICATION

Product Code	Priming Volume, mL	Flow Rate, mL/hr
A3060SL02	0.35	120
C3060SL5002	0.35	120
A4060SL02	0.62	540
C4060SL5002	0.62	540

PRODUCT CODE COMMENTS

A-30-60-SL-02

A >>> OTN Tearaway Introducer
30 >>> 3F
60 >>> Length
SL >>> Single Lumen
02 >>> NO Puncture Kit

C- 40-60-SL-50-02

C >>> Seldinger
40 >>> 4F
60 >>> Length
SL >>> Single Lumen
50 >>> Microintroducer Length
02 >>> NO Puncture Kit

A-OTN Tearaway Introducer

C-Seldinger

SL-Single Lumen


DL-Double Lumen

50-Microintroducer

02-NO Puncture Kit

AnderMed

Peripherally Inserted Central Venous Catheters [Distal Valve PU Material]

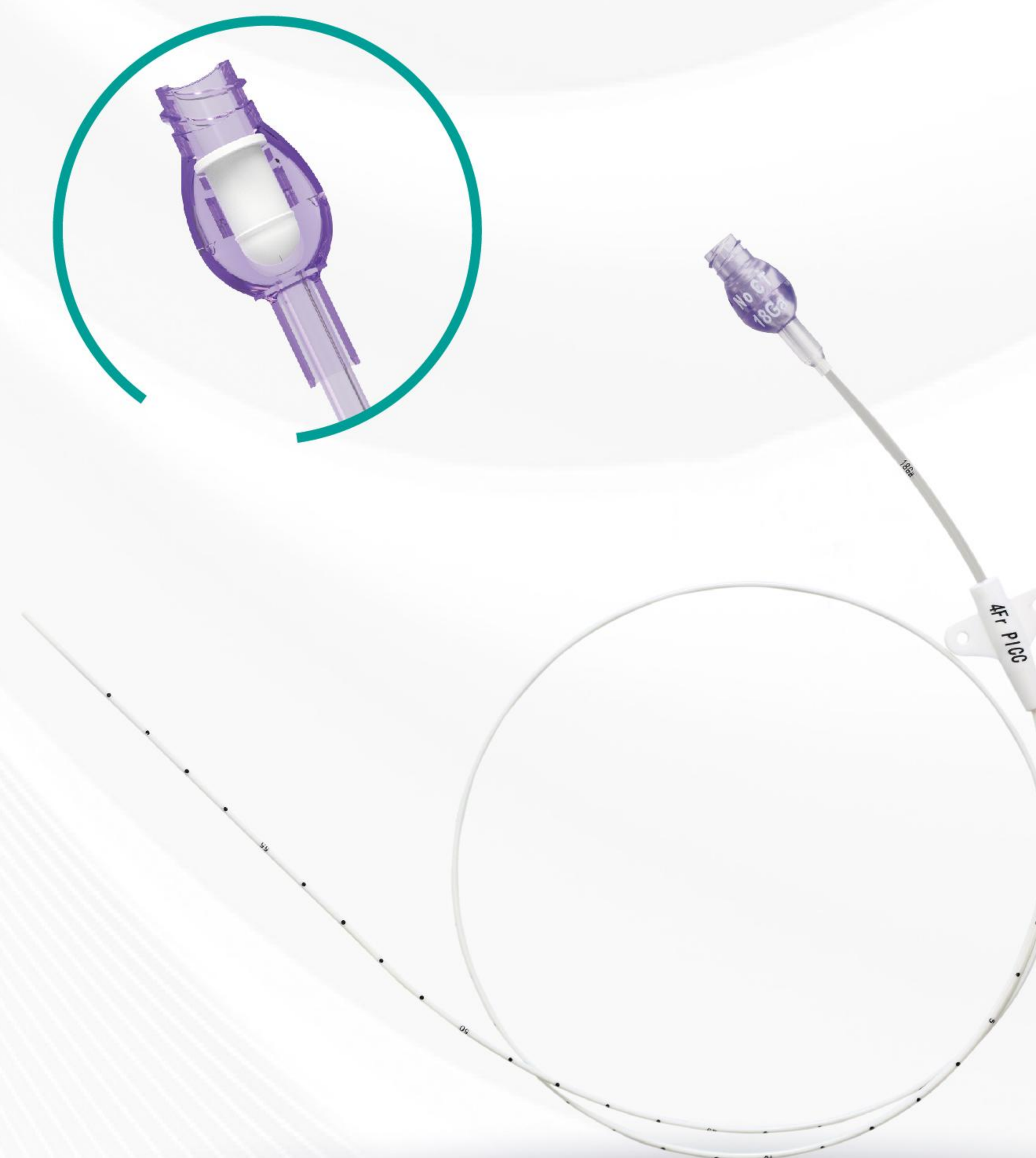

Shandong Ande Healthcare Apparatus Co., Ltd.



📍 Distal Valved Technology

- ⦿ Prevent blood reflux and reduce the occlusion of catheter, heparin is not needed.

- ✓ Valve opens allowing for infusion and flushing when positive pressure is applied.

- ✓ Valve opens allowing for aspiration when negative pressure is applied.

- ✓ The valve remains closed when it is not in use, reduce the risk of blood reflux and CRBSI.

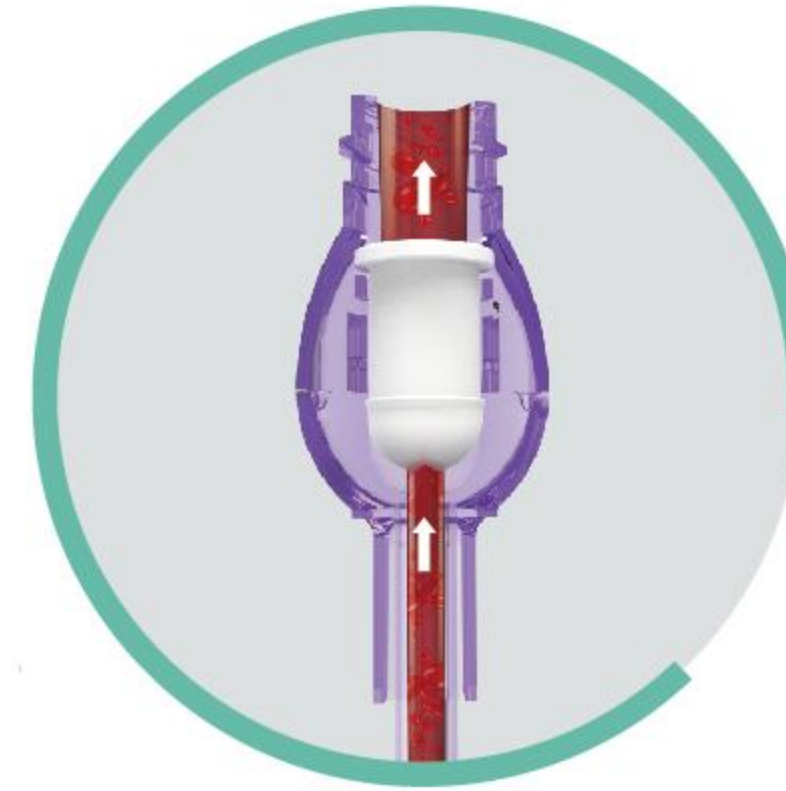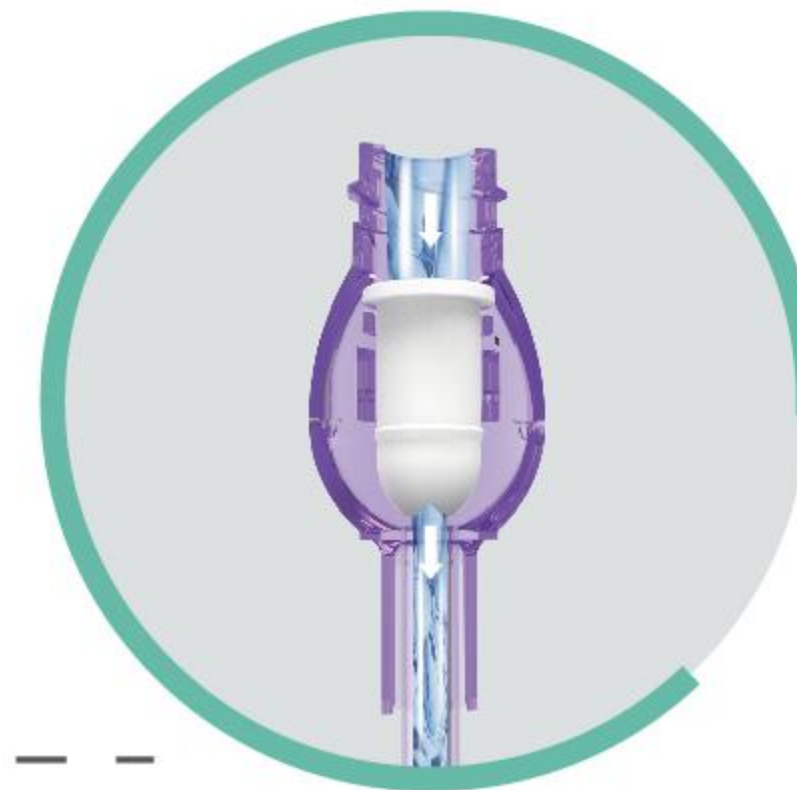

📍 Polyurethane Material

- ⦿ The catheter is flexible, tear and corrosion resistance, avoid leakage and breakage of catheter.
- ⦿ Smooth walls reduce adsorption, limit phlebitis, thrombosis and CRBSI.
- ⦿ Excellent biocompatibility, catheter softens with body temperature, better indwelling effect.



📍 Improved Seldinger kit

- ⦿ Improve the puncture success rate and reduce complications.

📍 Split-septum Neutral Needleless Connector

- ⦿ Reduce the risk of blood reflux and CRBSI.
Straight fluid path and clear housing enhances flushing efficiency and facilitates visualization of entire fluid channel.

📍 Integrated Design

- ⦿ Easy to use, avoid leakage and detachment of catheter.



4Fr PICC

SPECIFICATION

Product Code	Priming Volume, mL	Flow Rate, mL/hr
A3030SL02	0.3	192
C3030SL5002	0.3	192
A4030SL02	0.4	840
C4030SL5002	0.4	840
A4040SL02	0.5	840
C4040SL5002	0.5	840
A5040SL02	0.54	1380
C5040SL5002	0.54	1380
A5040DL02	0.46/0.46	600/600
C5040DL5002	0.46/0.46	600/600

PRODUCT CODE COMMENTS

A-30-60-SL-02

A >>> OTN Tearaway Introducer
30 >>> 3F
60 >>> Length
SL >>> Single Lumen
02 >>> NO Puncture Kit

C- 50-40-DL-50-02

C >>> Seldinger
50 >>> 5F
40 >>> Length
DL >>> Double Lumen
50 >>> Microintroducer Length
02 >>> NO Puncture Kit

A-OTN Tearaway Introducer

SL-Single Lumen

50-Microintroducer

C-Seldinger

DL-Double Lumen

02-NO Puncture Kit

AnderMed

Peripherally Inserted Central Venous Catheters [Midline Catheter]



Shandong Ande Healthcare Apparatus Co., Ltd.

📍 Distal Valved Technology

- 🎯 Prevent blood reflux and reduce the occlusion of catheter, heparin is not needed.

- ✓ Valve opens allowing for infusion and flushing when positive pressure is applied.

- ✓ Valve opens allowing for aspiration when negative pressure is applied.

- ✓ The valve remains closed when it is not in use, reduce the risk of blood reflux and CRBSI.

📍 Split-septum Needleless Connector

- 🎯 Reduce the risk of blood reflux and CRBSI.

Straight fluid path and clear housing enhances flushing efficiency and facilitates visualization of entire fluid channel.

📍 Multiple Lumens

- 🎯 High flow rate, limits infection rate, design for multiple clinical tasks: IV and blood administration, power injection, saline care and maintenance, etc.

📍 Integrated Design

- 🎯 Easy to use, avoid leakage and detachment of catheter.

📍 Power Injection Capability





- 🎯 Maximum injection rate of 5ml/s, maximum power injection pressure 300psi.
- 🎯 Universal catheter, design for power injection of contrast media and intravenous therapy.

📍 Improved Seldinger kit

- 🎯 Improve the puncture success rate and reduce complications.

📍 Polyurethane Material

- 🎯 The catheter is flexible, tear and corrosion resistance, avoid leakage and breakage of catheter.
- 🎯 Smooth walls reduce adsorption, limit phlebitis, thrombosis and CRBSI.
- 🎯 Excellent biocompatibility, catheter softens with body temperature, better indwelling effect.

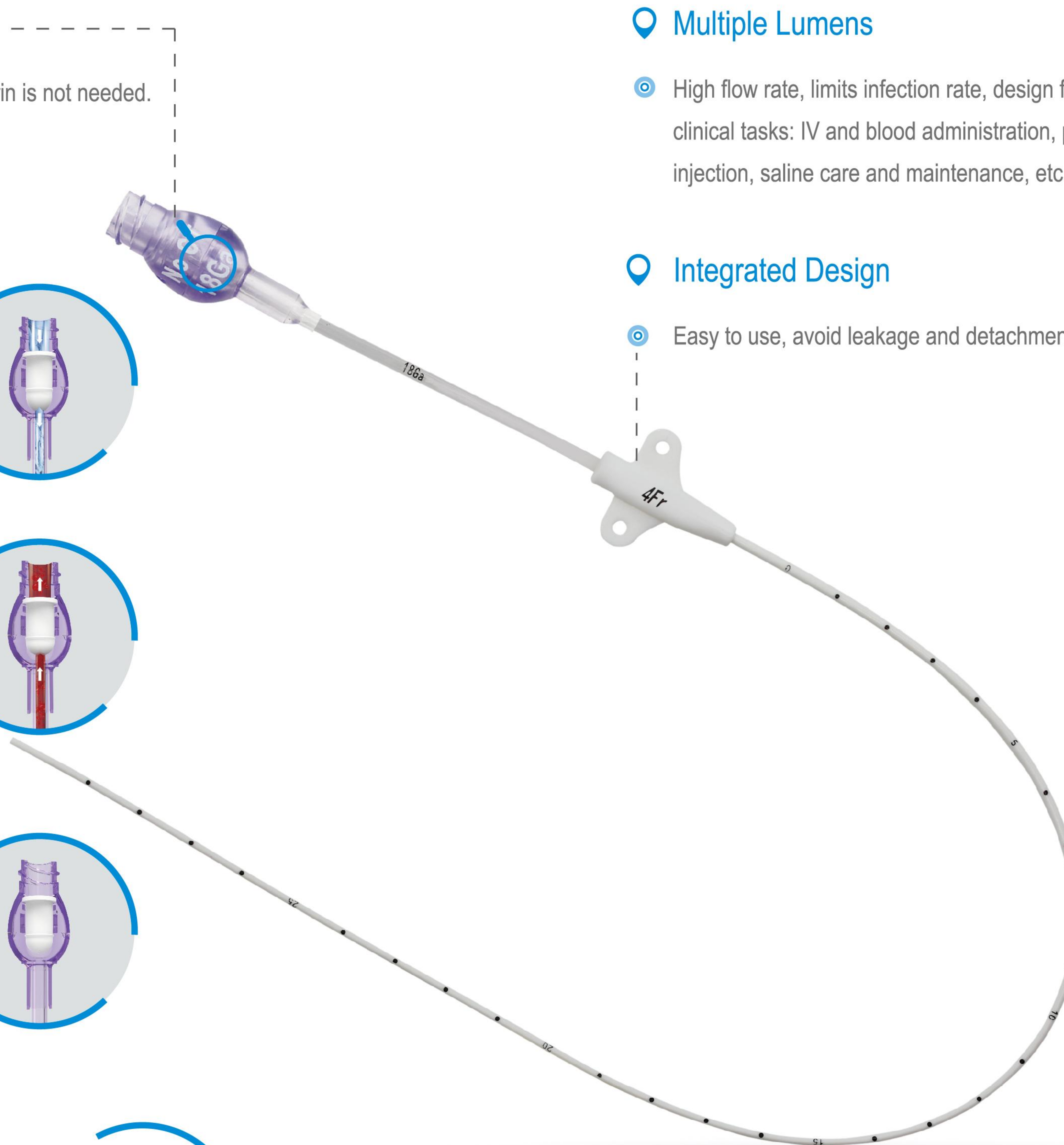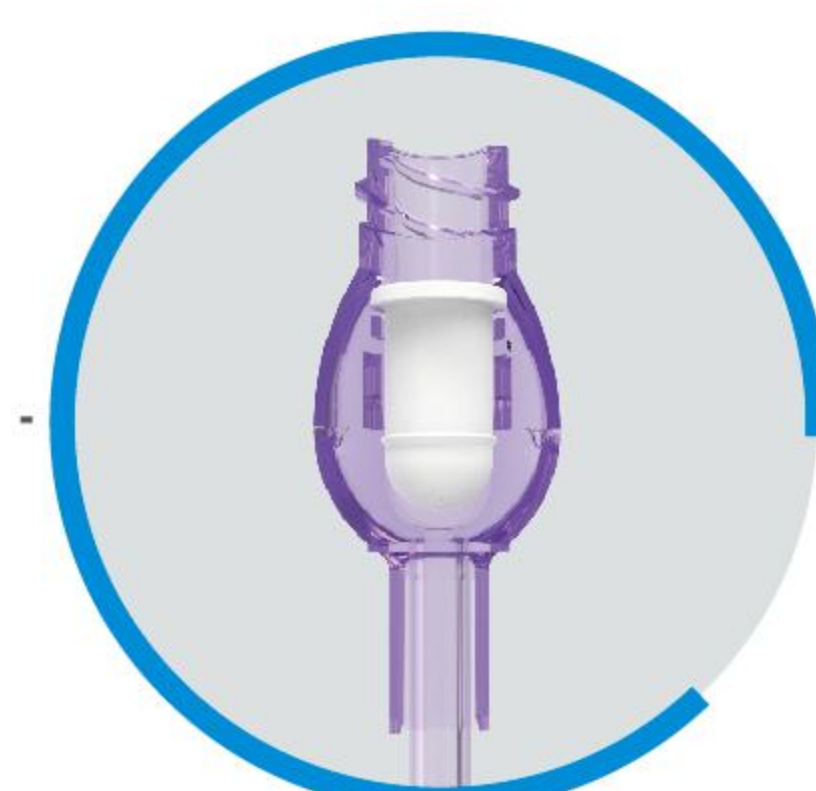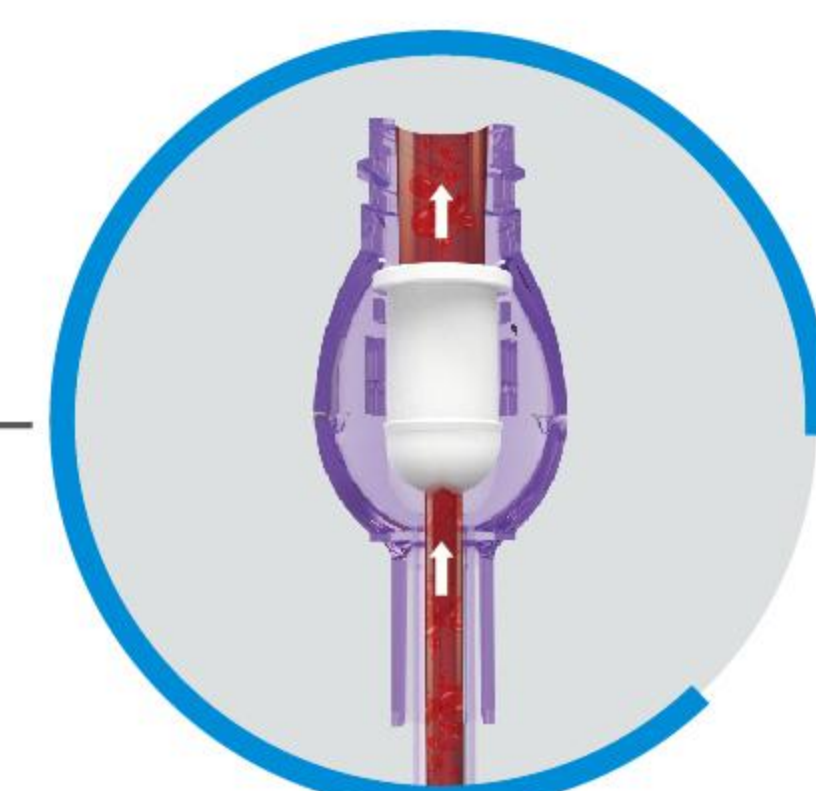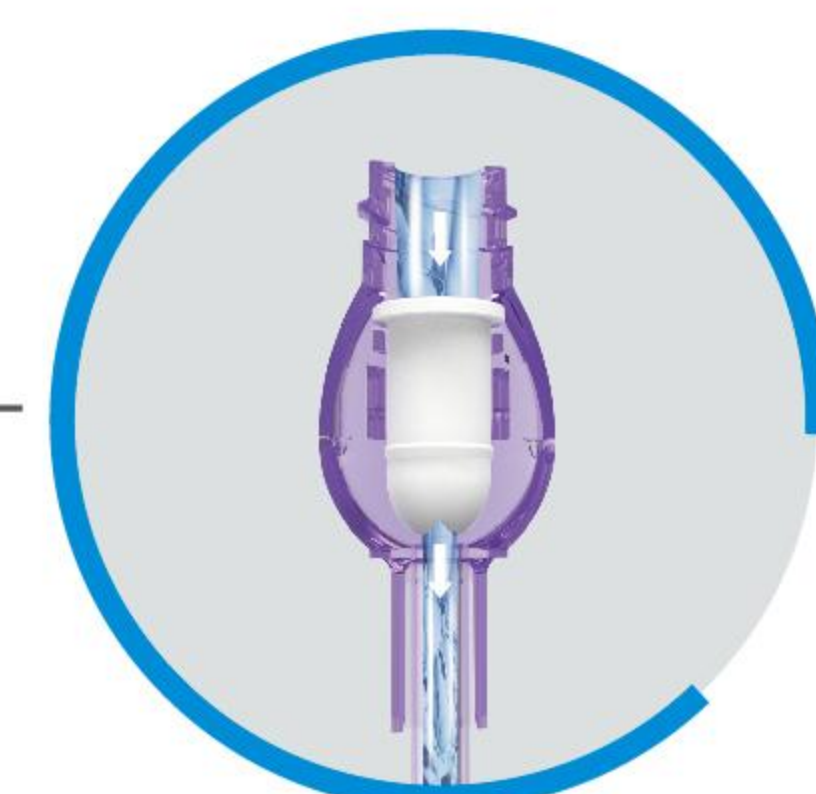

SPECIFICATION

Product Code	Priming Volume, mL	Flow Rate, mL/hr
A4065SL02	0.7	720
C4065SL5002	0.7	720
A5065SL02	0.75	930
C5065SL5002	0.75	930
A5065DL02	0.55/0.55	360/360
C5065DL5002	0.55/0.55	360/360
A4040SL02	0.5	840
C4040SL5002	0.5	840

Peripherally Inserted Central Venous Catheters
[CT USE]



PRODUCT CODE COMMENTS

C- 40-65-SL-50-02

C >>> Seldinger
40 >>> 4F
65 >>> Length
SL >>> Single Lumen
50 >>> Microintroducer Length
02 >>> NO Puncture Kit

C- 50-65-DL-50-02

C >>> Seldinger
50 >>> 5F
65 >>> Length
DL >>> Double Lumen
50 >>> Microintroducer Length
02 >>> NO Puncture Kit

A-OTN Tearaway Introducer

C-Seldinger

SL-Single Lumen

DL-Double Lumen

50-Microintroducer

02-NO Puncture Kit

Shandong Ande Healthcare Apparatus Co., Ltd.

Multiple Lumens


- High flow rate, limits infection rate, design for multiple clinical tasks: IV and blood administration, power injection, saline care and maintenance, etc.




Distal Valved Technology

- Prevent blood reflux and reduce the occlusion of catheter, heparin is not needed.


- ✓ Valve opens allowing for infusion and flushing when positive pressure is applied.

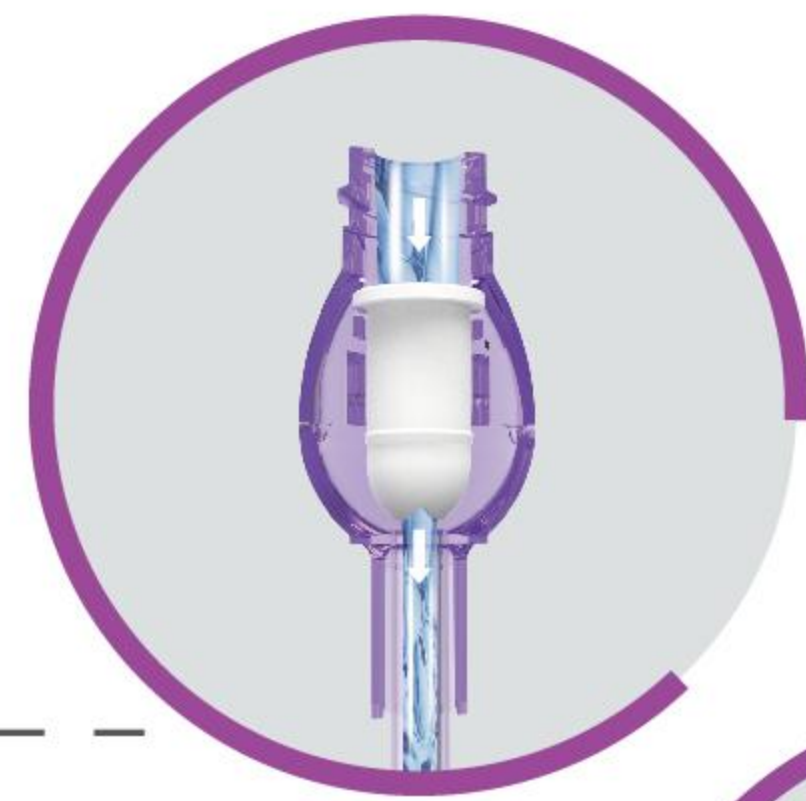

- ✓ Valve opens allowing for aspiration when negative pressure is applied.

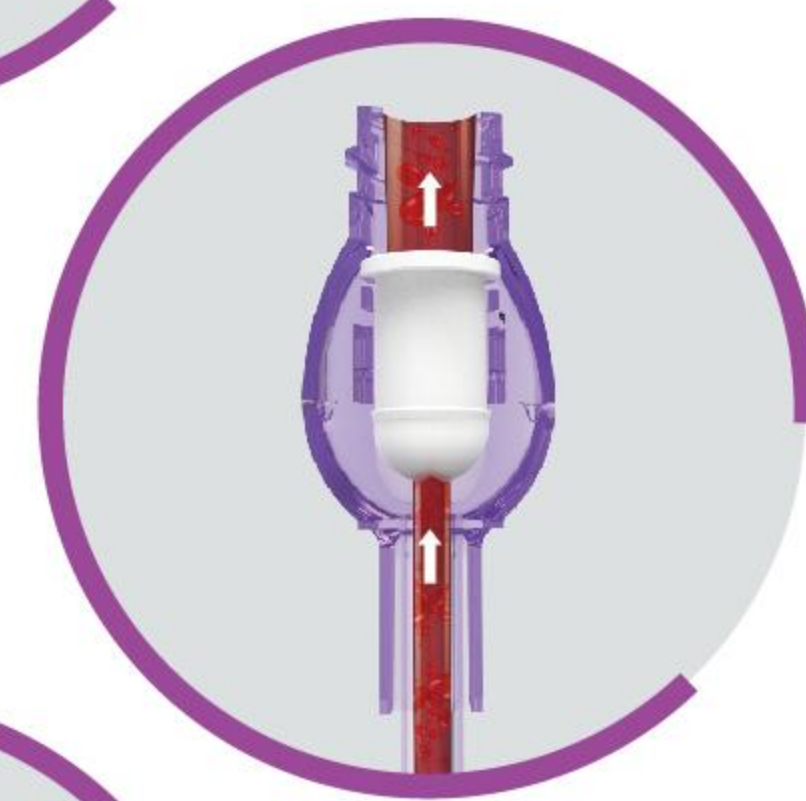

- ✓ The valve remains closed when it is not in use, reduce the risk of blood reflux and CRBSI.

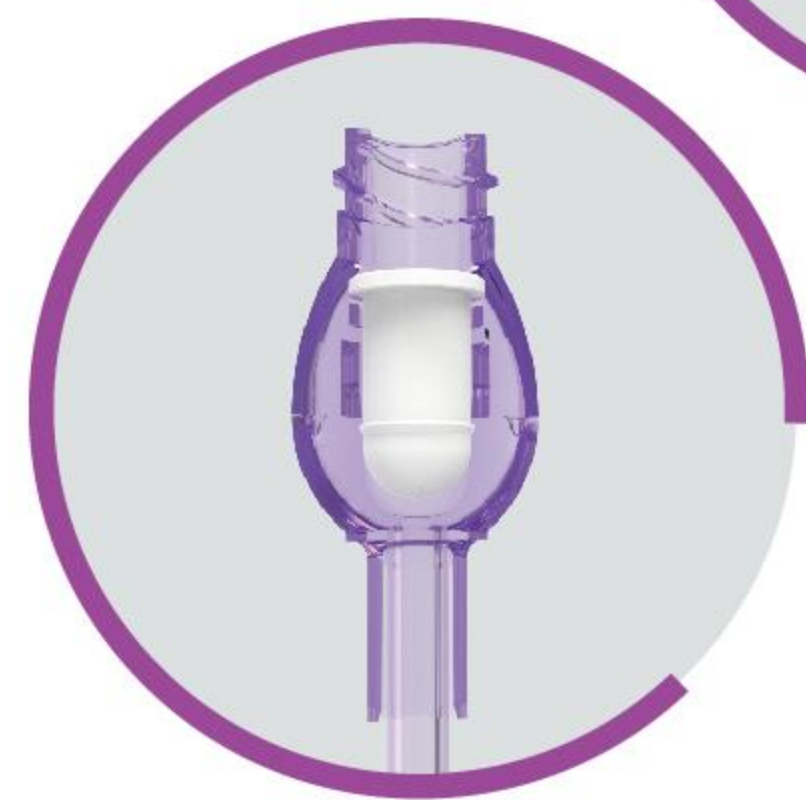

Split-septum Neutral Needleless Connector

- Reduce the risk of blood reflux and CRBSI. Straight fluid path and clear housing enhances flushing efficiency and facilitates visualization of entire fluid channel.



Integrated Design

- Easy to use, avoid leakage and detachment of catheter.

Reverse Tapered Catheter

- Fullfill the puncture site to avoid bleeding and lackage, easy to remove and kink resistant.

Polyurethane Material

- The catheter is flexible, tear and corrosion resistance, avoid leakage and breakage of catheter.
- Smooth walls reduce adsorption, limit phlebitis, thrombosis and CRBSI.
- Excellent biocompatibility, catheter softens with body temperature, better indwelling effect.



Improved Seldinger kit

- Improve the puncture success rate and reduce complications.

Power Injection Capability

- Maximum injection rate of 5ml/s, maximum power injection pressure 300psi.
- Universal catheter, design for power injection of contrast media and intravenous therapy.

SPECIFICATION

Product Code	Priming Volume, mL	Flow Rate, mL/hr
A1920SL02	0.17	84
A1950SL02	0.28	48

Peripherally Inserted Central Venous Catheters
[Neonate 1.9F]

PRODUCT CODE COMMENTS


A1920SL02	A1950SL02
A-OTN Tearaway Introducer	A-OTN Tearaway Introducer
19-1.9F	
20-Length	50-Length
SL-Single Lumen	
02-NO Puncture Kit	




📍 Distal Valved Technology

- 🕒 Prevent blood reflux and reduce the occlusion of catheter, heparin is not needed.


- ✓ Valve opens allowing for infusion and flushing when possitive pressure is applied. — — — —

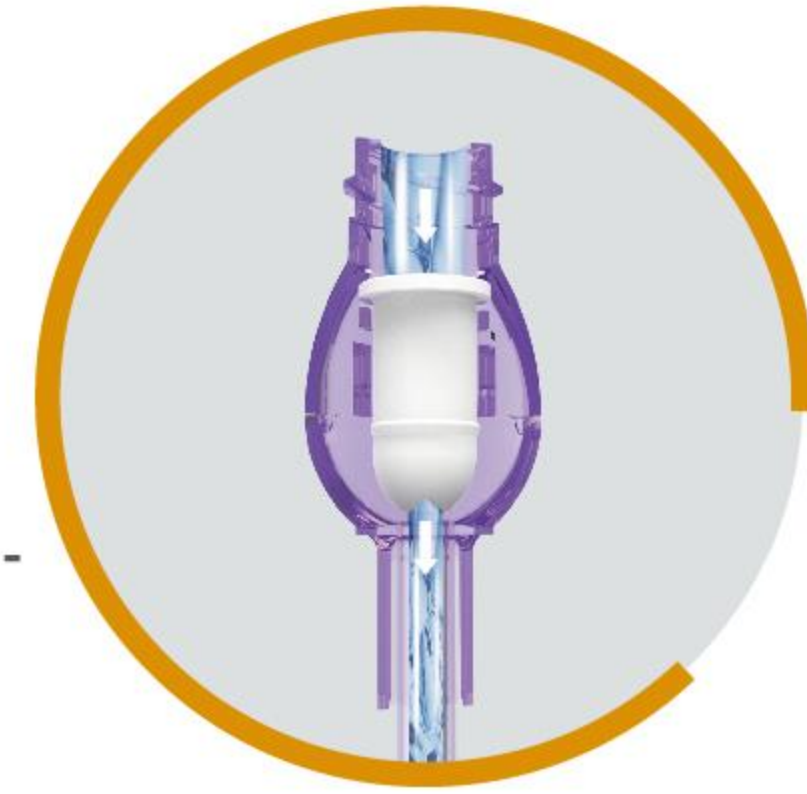

- ✓ Valve opens allowing for aspiration when negative pressure is applied. — — — —

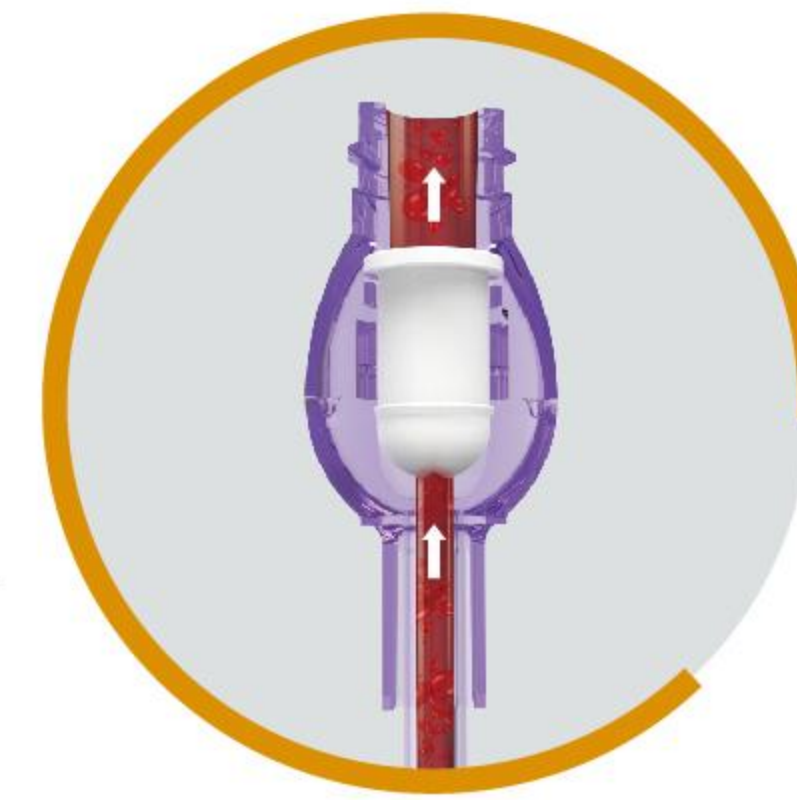

- ✓ The valve remains closed when it is not in use, reduce the risk of blood reflux and CRBSI. — — —

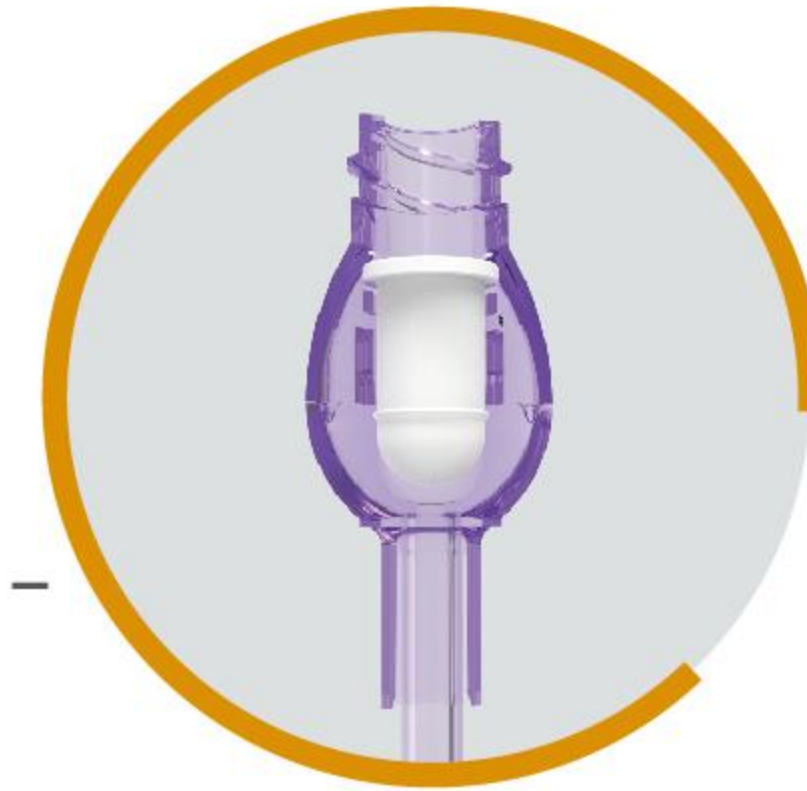

📍 Split-septum Neutral Needleless Connector

- 🕒 Reduce the risk of blood reflux and CRBSI.
Straight fluid path and clear housing enhances flushing efficiency and facilitates visualization of entire fluid channel.

📍 Integrated Design

- 🕒 Easy to use, avoid leakage and detachment of catheter.

📍 Polyurethane Material

- 🕒 The catheter is flexible, tear and corrosion resistance, avoid leakage and breakage of catheter.
- 🕒 Smooth walls reduce adsorption, limit phlebitis, thrombosis and CRBSI.
- 🕒 Excellent biocompatibility, catheter softens with body temperature, better indwelling effect.

